


SKIN+ACT

ARC LED Light Therapy Dome

PRODUCT MANUAL

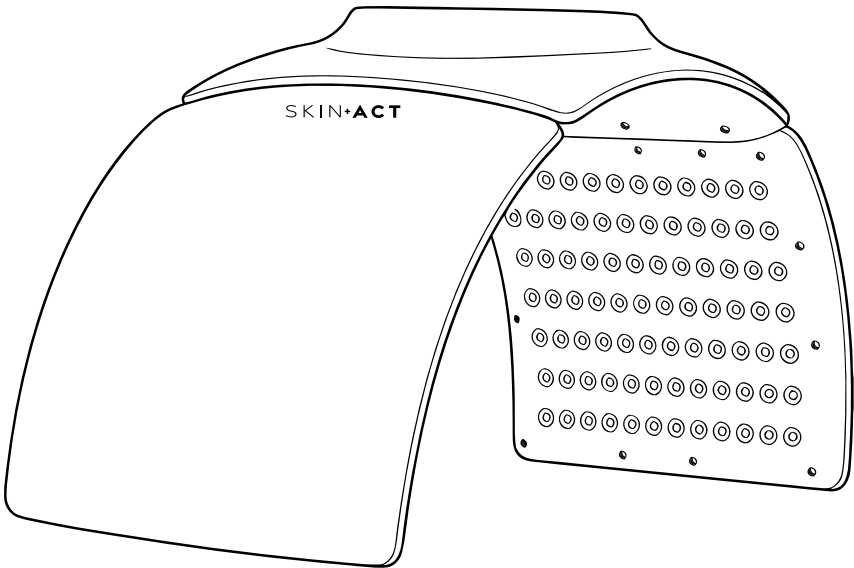

NOTICE

No part of this document may be photocopied, reproduced, or translated into any other language without prior written consent of Spa and Equipment/SkinAct.

The information in this document is subject to change without notice at any time and should not be construed as a commitment by Spa and Equipment/SkinAct.

Spa and Equipment/SkinAct assumes no responsibility for any errors that may be in this document, nor does it make any express or implied warranty with regard to this material, including, but not limited to, the implied warranties of merchantability and fitness for a particular purpose.

Spa and Equipment/SkinAct shall not be liable for incidental or consequential damages arising out of the furnishing, performance, or use of this document and the program material which it describes.

Thank You For Purchasing This Product!

Please read the entire instruction manual carefully before using the product. Take note of any warnings included in this manual.

Always keep this instruction manual handy and if you give the device to another person, hand over the manual as it is an essential part of the product.

Some components of the product may differ from the image that appears in catalogues, websites, brochures, etc.

In order to enhance the product, the manufacturer reserves the right to modify its features and/or components without prior notice.

Copyright © 2026 SkinAct, All Rights Reserved

TABLE OF CONTENTS

Disclaimer & Contraindications 4

Precautions & Safety Instructions 5

Included Items & Touch Panel Operation 6

Dual-Color Technology 7

Facial Protocol 8

Routine Maintenance & Intended Use 10

Specifications 11

DISCLAIMER

This product is intended for professional use by trained spa, salon, esthetic, wellness, or other qualified professionals. It must be used only for its intended purpose and in accordance with all instructions, warnings, contraindications, maintenance requirements, and safety guidelines provided in this manual.

The information in this manual is for general product operation and safety reference only. It does not replace professional training, licensing requirements, medical advice, service instructions, applicable laws, or workplace safety rules. Users are responsible for proper training, safe operation, inspection, client screening when applicable, sanitation, maintenance, and compliance with all applicable regulations.

SkinAct/Spa and Equipment does not provide medical advice and does not diagnose, treat, cure, or prevent any disease, medical condition, injury, or skin disorder. Any treatment descriptions, suggested uses, or expected results are for general cosmetic and professional use reference only, and results may vary.

Before each use, inspect the product, power cord, accessories, attachments, moving parts, and related components, as applicable, for damage, looseness, contamination, instability, or abnormal operation. Do not use the product if it appears damaged, defective, improperly assembled, unsafe, or if any component is missing or malfunctioning.

The user assumes full responsibility for proper assembly, installation, operation, client selection when applicable, treatment technique when applicable, cleaning, maintenance, and storage. SkinAct/Spa and Equipment shall not be responsible for injury, property damage, treatment complications, equipment damage, loss of business, or other damages resulting from misuse, improper use, lack of training, failure to follow instructions, unauthorized repairs or modifications, improper maintenance, or use outside the intended purpose.

Do not modify, disassemble, repair, or alter this product unless authorized in writing by SkinAct/Spa and Equipment or performed by an approved service provider. Unauthorized modification, misuse, improper maintenance, or failure to follow this manual may void the warranty.

Product designs, specifications, components, accessories, images, and instructions may vary from the actual product and may be changed or updated without prior notice.

Contraindications

DO NOT USE THIS DEVICE ON CLIENTS WHO:

- Are undergoing medical treatment or diagnostic procedures without authorization from a qualified healthcare professional.
- Have eye conditions or take photosensitizing medications without first consulting a qualified healthcare professional.

Precautions

- **Wear Eye Protection:** The operator and client must wear appropriate eye protection throughout the treatment.
- **Avoid Direct Eye Exposure:** Do not look directly into the LED light source while the device is operating.
- **Follow Treatment Times:** Typical treatment duration is 10–20 minutes, depending on the selected protocol and client tolerance.
- **Monitor Skin Response:** Mild warmth or temporary redness may occur. Discontinue treatment immediately if excessive heat, discomfort, irritation, or an adverse reaction develops.
- **Protect Skin After Treatment:** Apply sunscreen when the treated area may be exposed to sunlight.
- **Avoid Photosensitizing Products:** Do not use photosensitizing or skin-lightening ingredients during treatment unless they are approved for use with LED procedures.
- **Avoid Opaque Masks:** Opaque sheet masks may reduce or block light delivery and should not be used during treatment.
- **Inspect Before Use:** Check the device, attachments, and power cable for visible damage before each treatment.
- **Avoid Moisture:** Do not operate the device outdoors or in damp, wet, or excessively humid environments.


Safety Instructions

- Plug the device into a properly grounded electrical outlet with the appropriate voltage.
- Do not use the device if it or its power cable is damaged.
- Never touch the device, power plug, or cables with wet hands.
- Do not attempt to repair or modify the device. Repairs should only be performed by authorized technicians.
- Keep the device and accessories out of reach of children.

Included Items

Qty	Part Name	Image
1	ARC LED Light Therapy Dome	
1	Protective Eye Covers	
1	Power Adapter	

Touch Panel Operation

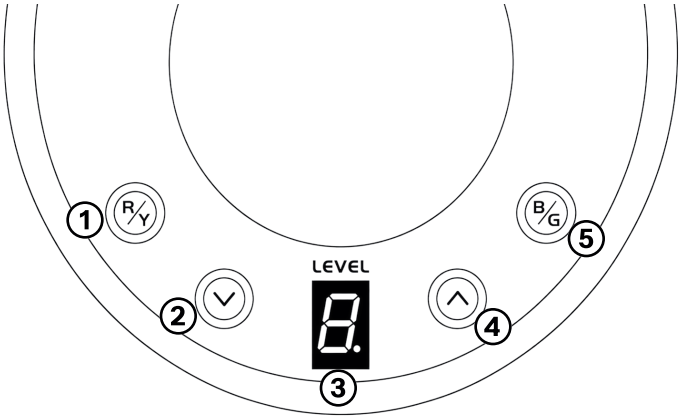

1. Red/Yellow Lights

2. Intensity Level Down

3. Digital Intensity Level Display
4. Intensity Level Up

5. Blue/Green Lights

Dual-Color Technology

The ARC Dome features a dual-side LED architecture:

- Side A: Blue or Green
- Side B: Red or Yellow

Available Combinations:

- **Red + Blue:** Supports collagen production and helps improve the appearance of blemish-prone skin.
- **Red + Green:** Supports collagen production and promotes a more even looking skin tone.
- **Yellow + Blue:** Helps soothe visible redness and improve the appearance of blemish-prone skin.
- **Yellow + Green:** Provides calming benefits and promotes a more balanced, even-looking complexion.

Each wavelength operates independently. Total energy delivered is additive.

Wavelength Characteristics

LED Color & Wavelength	Characteristics	Session Time
Blue (410 - 420 nm)	Lower-energy, surface-focused	10 - 20 min
Green (520 - 530 nm)	Moderate-energy, calming tone balance	10 - 20 min
Yellow (600-610 nm)	Moderate-energy, soothing protocols	10 - 20 min
Red (650-660 nm)	Higher-energy cosmetic interaction	10 - 20 min

Operating Instructions

1. Press R/Y to select either red or yellow light.
2. Press B/G to select either blue or green light.
3. For dual-color mode, select one color from each group to operate simultaneously.
4. Use the + / - controls to adjust the intensity.
5. Recommended treatment time: 10–20 minutes.

Light Therapy Facial Protocol

1. Cleanse

Thoroughly cleanse the face to remove makeup, sunscreen, oil, and surface debris. Perform a skin analysis and confirm that there are no contraindications present, ensuring LED treatment is appropriate for the client.

2. Exfoliate

Perform an appropriate gentle exfoliation based on client's skin type, condition, and treatment goals. Remove the exfoliant completely, using dampened facial sponges, and reassess the skin before proceeding.

3. Perform Extractions, if Needed

Complete gentle extractions when indicated. Do not perform extractions over irritated, broken skin or over any suspicious lesions. Remove any remaining product or debris and leave the skin clean before LED application.

4. Prepare for LED Therapy

Ensure the skin is clean and free of any opaque products that may block light delivery. The ARC LED Dome may be used on clean, bare skin or with a compatible, non-opaques facial serum applied beforehand.

Avoid photosensitizing or skin-lightening ingredients during LED exposure. Ensure the client is comfortably positioned and wearing the required protective eye covers. The licensed professional should also wear protective eyewear while the device is in operation.

Before placing the device over the client's face, **turn the ARC LED Dome on at its lowest intensity setting.** Once illuminated, carefully position the dome over the treatment area and then gradually adjust the intensity as appropriate for the client and treatment goals. Beginning at the lowest setting before positioning the device helps prevent the client from being startled by sudden activation of the LED lights, even while protective eyewear is in place.

5. Select LED Settings

Choose the appropriate LED color or dual-color combination according to the client's skin concerns and treatment goals. Select an appropriate intensity level, beginning conservatively when treating a new client.

Blue: 410–420 nm

Green: 520–532nm

Yellow: 600–660nm

Red: 650–660nm

6. Perform LED Treatment

Activate the ARC LED Dome and treat the skin for approximately **10–20 minutes**. Monitor the client throughout treatment and discontinue use if excessive warmth or discomfort occurs.

Professional Safety Note: The licensed professional should remain present and continue monitoring the client throughout the entire LED treatment. Staying nearby allows the professional to quickly respond if the client experiences discomfort, unexpected skin sensitivity, anxiety, eye-protection issues, or any other reaction that requires the treatment to be adjusted or stopped.

7. Final Product Application

Following LED therapy, apply the appropriate serum, followed by moisturizer and broad-spectrum SPF to hydrate, support, and protect the skin.

Treatment Time Limits

- **Typical Treatments:** 10–20 minutes
- **Facial Treatments:** Do not exceed 25 minutes
- **Body Treatments:** Do not exceed 30 minutes

Routine Maintenance

- Always unplug the unit before cleaning or performing maintenance.
- Wipe the exterior with a soft cloth lightly dampened with water and mild detergent.
- Do not immerse the unit in water or allow liquid to enter any opening.
- Allow the device to dry completely before use or storage.
- Store in a clean, dry, well-ventilated area.
- Have the device inspected annually by an authorized service technician.

Intended Use

This device is intended for professional cosmetic and esthetic use only. It is not intended to diagnose, treat, cure, or prevent any disease or medical condition. Operate the device in accordance with all applicable regulations and professional standards.

Specifications

Product Name:	ARC LED Light Therapy Dome
Voltage:	110-120 V AC
Frequency:	50/60 Hz
Power Consumption:	65 W
Adapter Output:	24 V DC, 3 A
Dimensions(L*W*H):	16 x 14 x 12.7”
LED Panels:	3 panels, each 10 × 10”
Intensity Levels:	1-9
Wavelengths:	Red, Blue, Yellow, Green
Total Number of LEDs:	432
Net Weight:	4.8 lb