

SKIN+ACT

Lucera Professional Electric Treatment/Procedure Chair (Table)

PRODUCT MANUAL

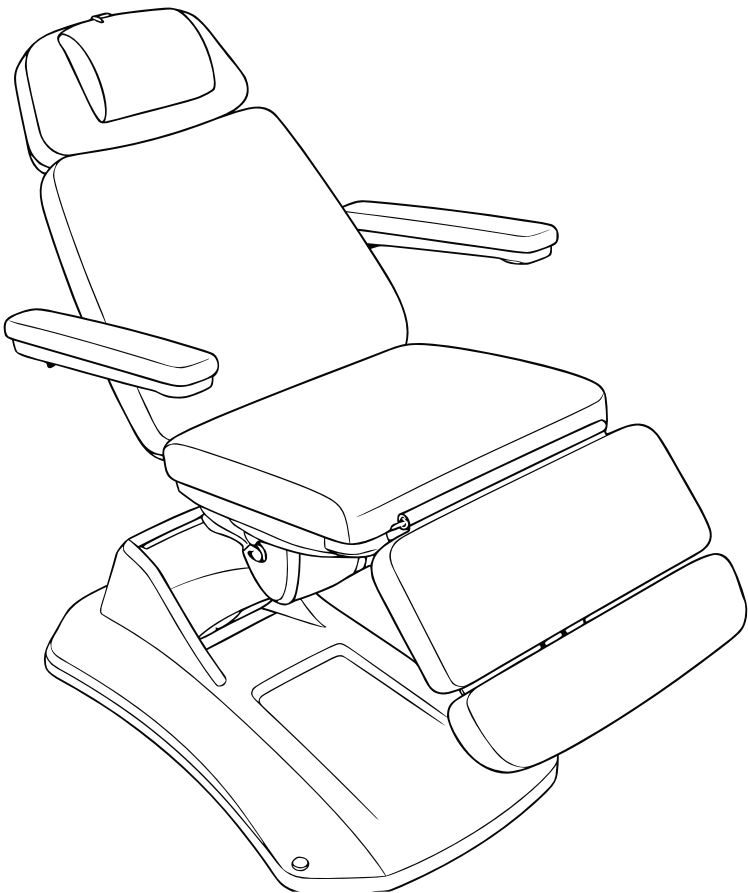


TABLE OF CONTENTS

Notice & Disclaimer3

Stability and Rotation Safety.....4

Safety and Operating Instructions.....5

Precautions.....6

Angle Specifications.....7

Dimensions & Structure of Treatment Table.....8

Installation.....9

Adjustment Controls on Backrest.....10

Adjustment Controls for Hand Controller.....11

Headrest Adjustments.....12

Lighting Adjustment.....13

Emergency Stop Switch.....13

Operation Cautions.....14

Maintenance.....14

Troubleshooting.....15

Specifications.....16

Warranty Information.....17

NOTICE

No part of this document may be photocopied, reproduced, or translated into any other language without prior written consent of Spa and Equipment/SkinAct.

The information in this document is subject to change without notice at any time and should not be construed as a commitment by Spa and Equipment/SkinAct.

Spa and Equipment/SkinAct assumes no responsibility for any errors that may be in this document, nor does it make any express or implied warranty with regard to this material, including, but not limited to, the implied warranties of merchantability and fitness for a particular purpose.

Spa and Equipment/SkinAct shall not be liable for incidental or consequential damages arising out of the furnishing, performance, or use of this document and the program material which it describes.

DISCLAIMER

This product is intended for professional use by trained spa, salon, esthetic, wellness, or other qualified professionals. It must be used only for its intended purpose and in accordance with all instructions, warnings, contraindications, maintenance requirements, and safety guidelines provided in this manual.

The information in this manual is for general product operation and safety reference only. It does not replace professional training, licensing requirements, medical advice, service instructions, applicable laws, or workplace safety rules. Users are responsible for proper training, safe operation, inspection, client screening when applicable, sanitation, maintenance, and compliance with all applicable regulations.

SkinAct/Spa and Equipment does not provide medical advice and does not diagnose, treat, cure, or prevent any disease, medical condition, injury, or skin disorder. Any treatment descriptions, suggested uses, or expected results are for general cosmetic and professional use reference only, and results may vary.

Before each use, inspect the product, power cord, accessories, attachments, moving parts, and related components, as applicable, for damage, looseness, contamination, instability, or abnormal operation. Do not use the product if it appears damaged, defective, improperly assembled, unsafe, or if any component is missing or malfunctioning.

The user assumes full responsibility for proper assembly, installation, operation, client selection when applicable, treatment technique when applicable, cleaning, maintenance, and storage. SkinAct/Spa and Equipment shall not be responsible for injury, property damage, treatment complications, equipment damage, loss of business, or other damages resulting from misuse, improper use, lack of training, failure to follow instructions, unauthorized repairs or modifications, improper maintenance, or use outside the intended purpose.

Do not modify, disassemble, repair, or alter this product unless authorized in writing by SkinAct/Spa and Equipment or performed by an approved service provider. Unauthorized modification, misuse, improper maintenance, or failure to follow this manual may void the warranty.

Product designs, specifications, components, accessories, images, and instructions may vary from the actual product and may be changed or updated without prior notice.

This product contains electric adjustment components and moving parts. Improper use may cause injury, pinching, equipment damage, or unsafe operation. Users are responsible for operating the table safely and following all instructions, warnings, and maintenance requirements in this manual. The user is solely responsible for confirming that use of this device is permitted under applicable state/local laws, scope-of-practice rules, licensing requirements, and workplace safety standards. This device is not intended for self-use, home use, or use by untrained individuals.

Safety and Operating Instructions

To ensure safe and proper operation of this electric treatment table, please read the following instructions carefully and retain this manual for future reference.

- Disconnect the power supply before replacing electrical components or performing cleaning or maintenance.
- Connect the power cord only to a properly installed and grounded outlet to reduce the risk of fire or electric shock.
- Do not stand, sit, or climb on any part of the equipment not intended to support the user.
- Keep children away from the equipment. Do not allow children to play with or operate the equipment.
- Before leaving the equipment unattended, return the backrest to its highest upright position.
- Do not install, operate, or store the equipment in wet or humid locations, such as bathrooms.
- Keep the hand controller, foot controller, and all electrical components away from water and other liquids.
- Do not operate the equipment near flammable anesthetics, oxygen-enriched environments, or other combustible gases.
- Individuals with physical, sensory, or cognitive limitations must use the equipment only under the supervision of a responsible person.
- Do not use the equipment if the power cord or plug is damaged. Contact the distributor, authorized service center, or a qualified technician for inspection and repair.
- The equipment contains no user-serviceable parts. Do not disassemble or attempt unauthorized repairs.
- Ensure there are no obstacles within the equipment's range of movement before making any adjustments.
- For indoor use only.

Stability and Rotation Safety

When operating a hydraulic rotating chair or single-post electric treatment table, the operator is responsible for maintaining stability and ensuring client safety. Always follow the instructions provided for the specific model.

To ensure stability

1. Place the chair or table on a level, dry, and non-slip surface. Adjust the leveling feet, if equipped, so the base remains stable and does not rock.
2. Inspect the base, controls, adjustment mechanisms, and fasteners before each use. Do not operate the unit if any component is loose, damaged, or functioning improperly.
3. Keep the client's weight centered over the base. Position the chair or tabletop in line with the base whenever possible, particularly during treatment.
4. Do not allow anyone to sit, stand, lean, or place excessive weight on the headrest, backrest, leg rest, foot section, or armrests.
5. Whenever possible, rotate or reposition the unit before the client enters it. Do not rotate or adjust the unit while the client is entering or exiting.
6. If rotation is necessary while occupied, place the client in a stable upright position, keep the client centered, support the chair or table securely, and rotate it slowly.
7. Keep hands, feet, clothing, cords, and other objects away from moving parts and potential pinch points during rotation or adjustment.
8. Ensure all rotation locks and adjustable sections are securely engaged before beginning treatment.
9. Do not exceed the stated weight capacity or apply sudden, concentrated, or excessive downward pressure. These products are not intended for massage or procedures that require substantial downward force.
10. Do not leave a client unattended while the unit is elevated, reclined, rotating, or being adjusted.
11. Stop using the unit immediately if it becomes unstable, shifts unexpectedly, makes unusual sounds, or fails to hold its selected position.



Warning: Do Not Sit or Apply Excess Force



Do not sit, stand, jump, or place excessive weight on the backrest or footrest. Doing so may cause the product to tip, damage the backrest/footrest joints, or result in injury.

30%

Back Cushion

55%

Seat Cushion

15%

Leg Cushion






Precautions

- This product must be operated by trained and experienced personnel who are familiar with its structure, functions, and operating instructions.
- Read this manual carefully before installation, operation, or adjustment. Keep the manual available for future reference.
- Do not stand, sit, climb, or place excessive weight on any part of the product not designed to support a person.
- Follow all operating limits and safety instructions. Do not overload the product or modify its safety features.
- Inspect the product regularly and tighten any loose bolts, nuts, or fasteners to prevent movement, vibration, or component failure.
- Keep hands, feet, hair, clothing, and other body parts away from moving components during adjustment.
- Do not operate the product if it is damaged, malfunctioning, or producing unusual sounds. Stop use immediately and contact an authorized service technician.
- Operate and store the product only in a clean, dry environment that meets the specified operating requirements.
- Connect the three-prong power cord only to a properly installed and grounded three-prong outlet. Do not remove or bypass the grounding pin.
- Follow all warning labels attached to the product.
- Press and hold the appropriate control button while adjusting the product. Release the button when the desired position is reached.
- To prevent the motor from overheating, do not operate it continuously for more than 2 minutes. Allow the motor to rest for at least 18 minutes before operating it again.
- Maintenance and repairs must be performed only by trained professionals.
- After each treatment, clean the product with a mild, neutral cleaning solution and a soft cloth.
- Keep hands, feet, and head clear of all pinch points and moving areas while adjusting the product.


Angle Specifications

Rotation
-120° to +120°

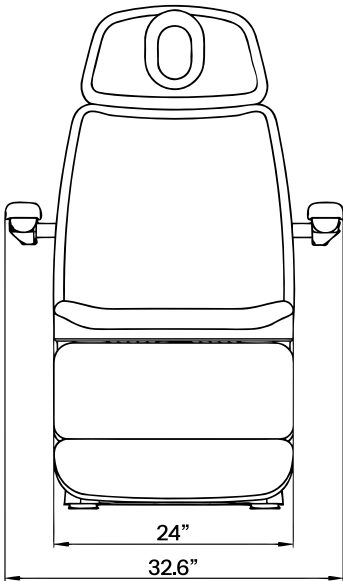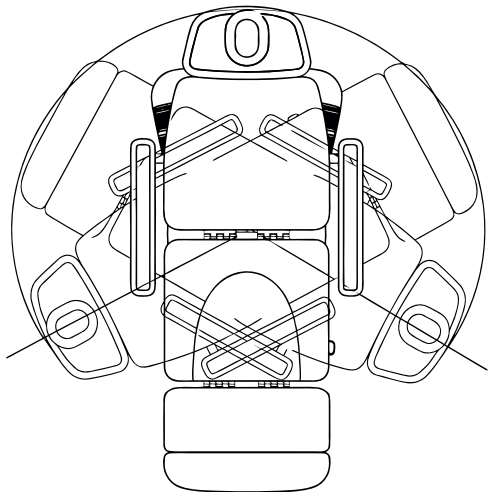

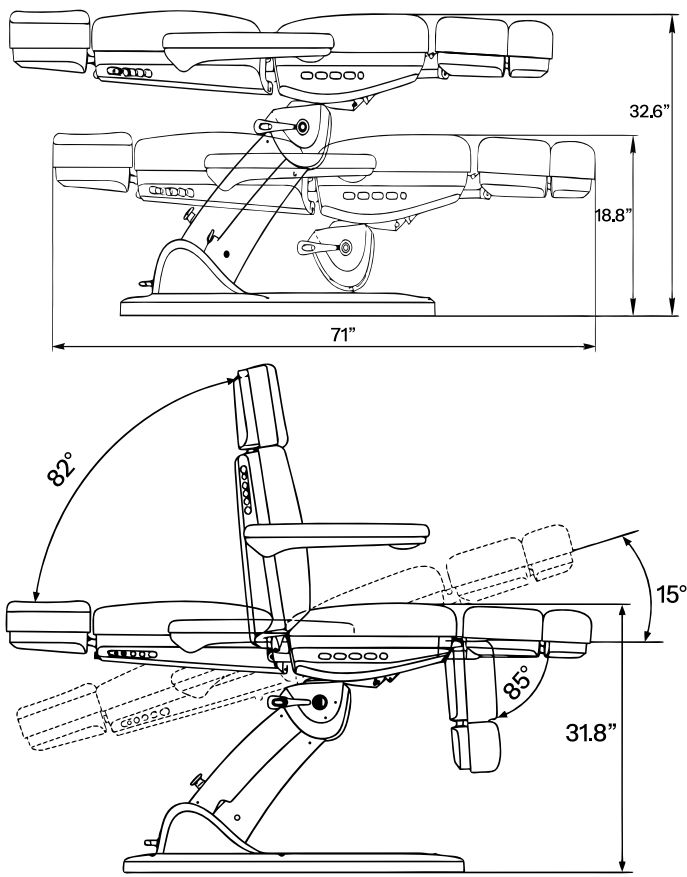

Structure of Treatment Table

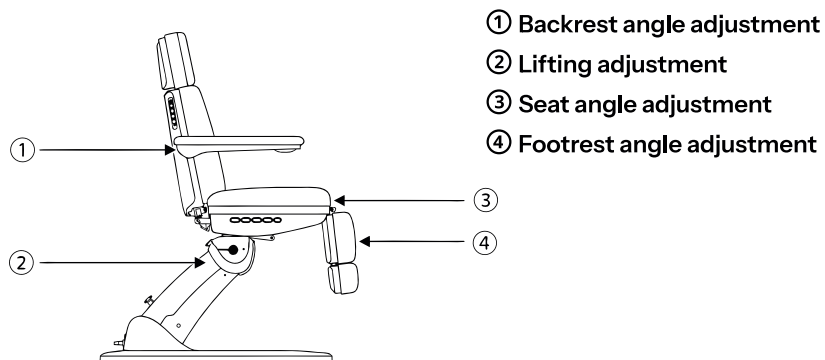

Installation

1. Choose an appropriate location.

Place the treatment table in a clean, dry, and cool environment. Select a well-lit location with enough space for safe operation and easy access around the table.

2. Inspect the product before installation.

Open the packaging and confirm that the treatment table is undamaged. Check that all accessories, spare parts, and components are included. Contact the distributor or SkinAct if anything is missing or damaged.

3. Position the treatment table.


Place the table on a firm, flat, and stable surface. No additional floor installation or fastening is required, as the table is designed to remain stable when positioned correctly.

4. Level the table if needed.


If the table rocks or appears uneven, adjust the four plastic leveling feet located underneath the base. Use an open-end wrench to adjust the feet and M10 locking nuts until the table is level and stable.

5. Connect the power plug.

Connect the power plug securely to a properly grounded power outlet before operating the table.

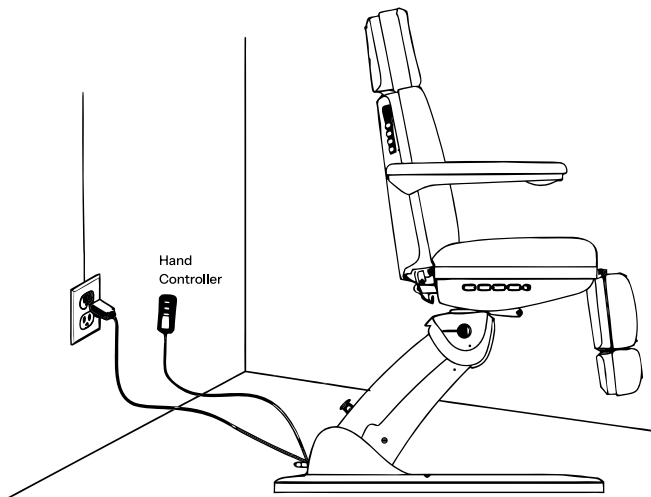


Adjustment Controls on Backrest

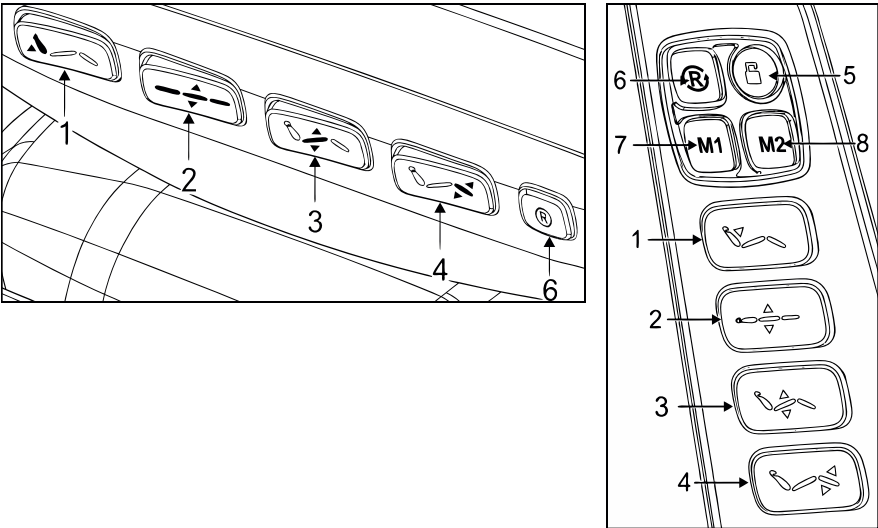

No.	Adjustment	Operation
1	Backrest Adjustment	Press and hold the Up or Down button to raise or lower the backrest. The backrest angle adjusts at a constant speed.
2	Height Adjustment	Press and hold the Up or Down button to raise or lower the entire treatment table. The table moves at a constant speed.
3	Seat Tilt Adjustment	Press and hold the appropriate button to tilt the seat cushion forward or backward. The seat angle adjusts at a constant speed.
4	Leg Rest Adjustment	Press and hold the Up or Down button to raise or lower the leg rest. The leg rest angle adjusts at a constant speed.
5	Safety Lock Button	Press and hold the Safety Lock button for approximately 3 seconds until the red indicator light turns on. All control buttons will be disabled except the Safety Lock button. Press and hold the button again for approximately 3 seconds to unlock the controls and return them to normal operation.
6	One-Touch Reset Button	Press the button to automatically lower the treatment table to its lowest position.
7/8	Memory Position 1 & 2 – M1 & M2	Adjust the treatment table to the desired position. Press and hold the M1 or M2 button for 2–3 seconds. When the table beeps, the position has been successfully saved.


Adjustment Controls for Hand Controller

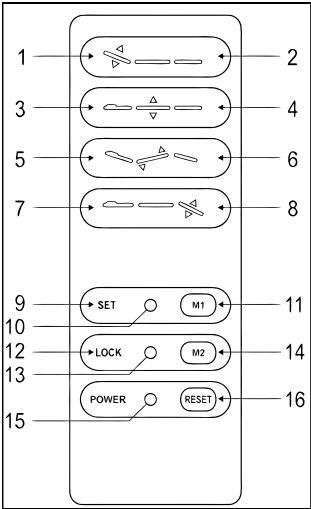


No.	Adjustment	Operation
1	Backrest Up	Press and hold the button to raise the backrest smoothly.
2	Backrest Down	Press and hold the button to lower the backrest smoothly.
3	Table Up	Press and hold the button to raise the entire treatment table smoothly.
4	Table Down	Press and hold the button to lower the entire treatment table smoothly.
5	Seat Cushion Up	Press and hold the button to tilt the seat cushion upward.
6	Seat Cushion Down	Press and hold the button to tilt the seat cushion downward.
7	Leg Rest Up	Press and hold the button to raise the leg rest smoothly.
8	Leg Rest Down	Press and hold the button to lower the leg rest smoothly.
9	Memory Setting Button	Press and hold the SET button for approximately 2 seconds to enter memory-setting mode.
10	Memory Setting Indicator	The indicator light turns on when memory-setting mode is activated.
11	Memory Position 1—M1	Press and hold SET for approximately 2 seconds. When the indicator light turns on, adjust the table to the desired position and press M1 to save it. Press M1 during normal operation to move the table automatically to the saved position.

No.	Adjustment	Operation
12	Safety Lock Button	Press and hold the Safety Lock button for approximately 2 seconds until the red indicator light turns on. All buttons will be disabled except the Safety Lock button. Press and hold it again for approximately 2 seconds to unlock the controls.
13	Safety Lock Indicator	The red indicator light turns on when the controls are locked and turns off when normal operation is restored.
14	Memory Position 2—M2	Press and hold SET for approximately 2 seconds. When the indicator light turns on, adjust the table to the desired position and press M2 to save it. Press M2 during normal operation to move the table automatically to the saved position.
15	Power Indicator	The indicator light turns on when the hand controller is properly connected and receiving power. If the light does not turn on, check the power connection.
16	One-Touch Reset Button	Press the button to automatically return the treatment table to its lowest position.

Headrest Adjustments

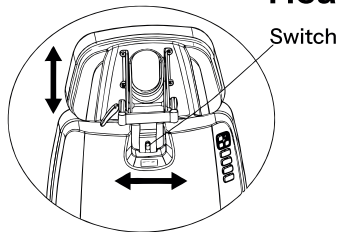


Turn the switch to the right to unlock the headrest. Adjust the headrest to the desired height, then turn the switch to the left to lock it securely in place.

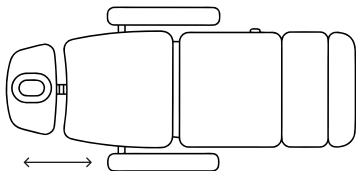


1. Shift the headrest up about 2".

2. Remove the backrest cushion as shown in the image.

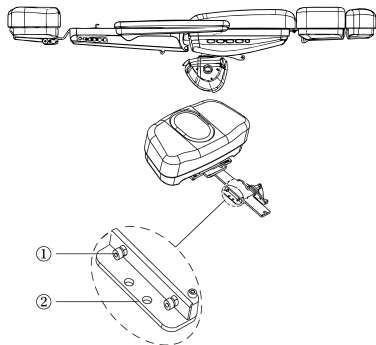


3. Use an open-end wrench and a hex wrench to tighten and secure the screws in positions ① and ②, as shown.

Lighting Adjustment

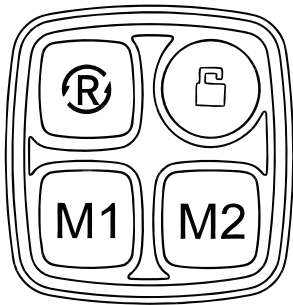

1. Briefly press the Lighting button to enter lighting-adjustment mode. The button indicators will flash rapidly.
2. Press **M1** to cycle through the **16 available lighting colors and modes**.
3. Press **M2** to adjust the brightness through **12 levels**.
4. After selecting the desired color and brightness, press the **Lighting button** again to exit lighting-adjustment mode.


If no button is pressed for 10 seconds, the controller will automatically exit lighting-adjustment mode.

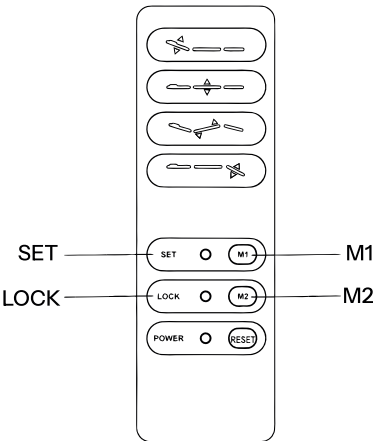

1. Press and hold the **SET/LOCK** button for approximately **2 seconds** to enter lighting-adjustment mode.
2. Press **M1** to cycle through the **16 available lighting colors and modes**.
3. Press and hold **M2** to adjust the brightness through **12 levels**.
4. After selecting the desired lighting mode and brightness, press and hold **SET/LOCK** to save the settings and exit lighting-adjustment mode.


If no button is pressed for 10 seconds, the controller will automatically exit lighting-adjustment mode.

Emergency Stop Switch

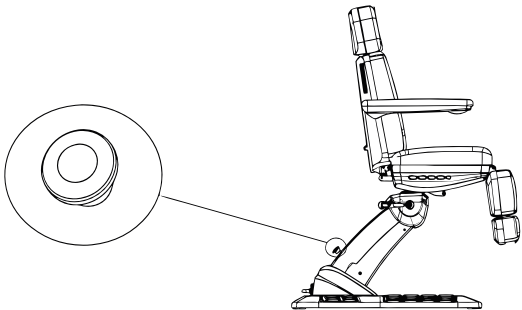

Immediately stops all treatment table movement in the event of an emergency or unexpected situation.

Operation Cautions

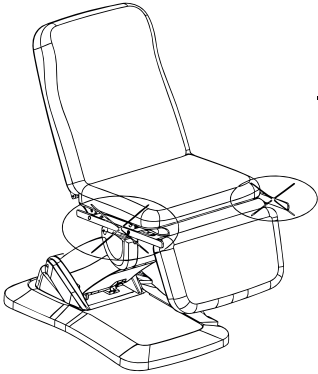

- Do not use the sliding guide as a handle or lifting point, as this may damage or break the component.

▪ This product graphic is for reference only. Please refer to the actual product for the specific structure.



WARNING
Please start the reset function after completely closing the foot extension.



NO SITTING



Please do not sit on the backrest/footrest or jump to the backrest/footrest, avoid the risk of falling due to excessive force fracture of the backrest/footrest joint

Maintenance

- When cleaning and sanitizing, do not use bleach, Barricade, or other harsh chemicals. Using harsh chemicals can damage the upholstery.
- For best results, we recommend cleaning with **SkinAct All-In-One Upholstery & Disinfecting Multi-Purpose Cleaner** (SKU24261).
- Clean and sanitize the upholstery after each treatment.
- Every six months, visually inspect all visible parts of the unit. If there are any abnormal noises or vibrations during operation, call Spa & Equipment/SkinAct at 1-800-418-0830.

Troubleshooting

Motor Does Not Operate Normally

Solution: Confirm that the power cord is securely connected to the treatment table and a properly grounded outlet. Check that the hand controller or foot controller is connected correctly. Make sure the emergency stop switch has been released and reset. If the motor still does not operate, the controller, control box, motor cable, or motor may be damaged. Stop using the table and contact SkinAct/Spa & Equipment for assistance.

Abnormal Noise During Operation

Solution: Inspect the hinges, lifting mechanism, and other moving components for dryness, looseness, or visible damage. Apply an appropriate lubricant to the designated moving points when needed. If the noise continues or appears to come from the motor, stop using the table and contact SkinAct/Spa & Equipment for service.

If the problem remains after completing the recommended checks, discontinue use and contact SkinAct/Spa & Equipment for further assistance.

Specifications

Product Name:	Lucera Electric Medical Spa Treatment Table/Chair
Voltage:	AC 100-240 V
Frequency:	50/60 Hz
Power Consumption:	300 VA max
Backrest Tilt Angle:	0° flat to 82° incline
Height Range:	18.8” - 32.6”
Length with Headrest:	71.5”
Table Tilt Angle Back:	0° flat to 15°
Leg Rest Angle Down:	0° flat to 85°
Width with Armrest:	32.6”
Width without Armrest:	24”
Weight Capacity:	440 lbs

Warranty Information

The manufacturer/warrantor warrants to the original purchaser that this product will be free from defects in materials and workmanship under normal and proper use for a period of one (1) year from the original date of purchase, unless otherwise stated in writing for a specific product.

This limited warranty applies only to the original purchaser and is not transferable. Proof of purchase may be required to obtain warranty service.

This warranty covers parts and labor required to correct defects in materials or workmanship during the warranty period. At the warrantor’s discretion, the product may be repaired, replaced, or substituted with a comparable item or component.

This warranty does not cover damage or defects caused by misuse, abuse, improper assembly, improper installation, unauthorized modification, unauthorized repair, neglect, accident, improper cleaning, use of harsh chemicals, power surges, water or liquid damage, use outside the intended purpose, failure to follow the product manual, or use with parts/accessories not approved by SkinAct.

This warranty does not cover normal wear and tear, cosmetic damage, upholstery wear, stains, discoloration, scratches, dents, wheels/casters, bulbs, glass parts, filters, tips, probes, disposable parts, consumables, or accessories subject to regular replacement unless the defect is determined to be due to a manufacturing issue.

The warrantor shall not be liable for loss of use, loss of time, loss of business, lost revenue, inconvenience, or any indirect, incidental, special, or consequential damages arising from the use of this product or from any defect covered by this warranty.

To request warranty service, contact the dealer or SkinAct customer support. If return service is required, the purchaser may be asked to return the product or the affected component for inspection. The purchaser is responsible for properly packaging the item and prepaying all shipping costs to send the item in for warranty service. If the issue is covered under warranty, repair or replacement will be performed at no charge, subject to the terms of this limited warranty.

Warranty does not cover damage, contamination, upholstery damage, sanitation-related issues, injury claims, or equipment issues resulting from improper cleaning, improper maintenance, misuse, use by untrained users, or use outside applicable safety guidelines. This limited warranty is the only written or express warranty provided by the warrantor. No employee, representative, dealer, or third party is authorized to modify or extend this warranty unless confirmed in writing by SkinAct. This warranty gives you specific legal rights, and you may also have other rights, which vary by state or province.

SKIN+ACT

Phone:1-800-418-0830

Email:info@skinact.com

Website:www.skinact.com