
SKIN+ACT

2 in 1 High Frequency & Galvanic Elite Edition

PRODUCT MANUAL

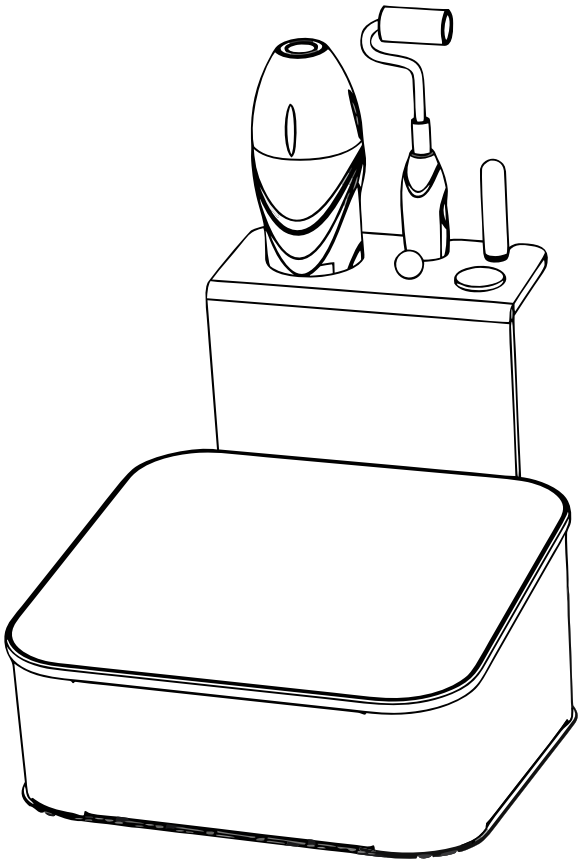

SKU25344

SKINACT.COM

NOTICE

No part of this document may be photocopied, reproduced, or translated into any other language without prior written consent of Spa and Equipment/SkinAct.

The information in this document is subject to change without notice at any time and should not be construed as a commitment by Spa and Equipment/SkinAct

Spa and Equipment/SkinAct assumes no responsibility for any errors that may be in this document, nor does it make expressed or implied warranty with regard to this material, including, but not limited to, the implied warranties of merchantability and fitness for a particular purpose.

Spa and Equipment/SkinAct shall not be liable for incidental or consequential damages arising out of the furnishing, performance, or use of this document and the program material which it describes.

Thank You For Purchasing This Product!

Please read the entire instruction manual carefully before using the product. Take note of any warnings included in this manual.

Always keep this instruction manual handy and if you give the device to another person, hand over the manual as it is an essential part of the product.

Some components of the product may differ from the image that appears in catalogues, websites, brochures, etc..

In order to enhance the product, the manufacturer reserves the right to modify its features and/or components without prior notice.

Copyright © 2026 SkinAct, All Rights Reserved

TABLE OF CONTENTS

Disclaimer & Contraindications 4

Contraindications & Safety Instructions 5

Warnings (High Frequency & Galvanic Mode) 6

Included Items 7

High Frequency Tool Components 8

Galvanic Tool Components 9

Touch Panel Interface 10

Device Components 11

Assembly Instructions 12

Operating Instructions 13

High Frequency Information 14

High Frequency Operating Instructions 15

Galvanic Operating Instructions 18

Usage (Iontophoresis Function) 19

Specifications 21

Warranty Information 22

DISCLAIMER

- This product has been designed to be used by esthetic professionals at spas and salons, it should not be used for other purposes.
- Please note that any misuse of the device entails the loss of the guarantee.
- Do not make any modification to the device without the authorization.
- We shall not be held liable for any material or personal damage due to the misuse of the product.
- Documents are made available without warranty of any kind. They may be updated, altered, expanded, or removed at any time. Skin Act does not warrant or guarantee any of the processes described herein or perform any independent analysis in connection with guidelines provided.
- Readers are expressly warned to consider and adopt all safety precautions indicated herein and to avoid all potential hazards. By following the instructions provided, the reader willingly assumes all risks in connection with such instructions.
- If the reader has questions about any of the guidelines provided by SkinAct we recommend they contact their nearest Esthetician Training Center to obtain more explicit instruction on Facial Machine Usage.

Contraindications

THIS DEVICE SHOULD NOT BE USED ON CLIENTS WITH:

- Pacemakers, metallic implants, or other electronic medical devices.
- Heart conditions or a history of heart disease.
- Pregnancy.
- Inflamed, irritated, or damaged skin.
- Severe acne rosacea or vascular conditions.
- Epilepsy or other seizure disorders.
- Clients undergoing medical treatments or with underlying health conditions must obtain prior clearance from a physician.

Precautions

- **Inspect the Device:** Check for any visible damage to the device or its power cable before and after each use.
- **Unplug When Not in Use:** Always disconnect the device from the power source before cleaning, storing, or performing maintenance.
- **Avoid Humidity:** Do not operate the device in bathrooms or other excessively humid environments.
- **Position Cables Safely:** Ensure the power cable is positioned to prevent tripping hazards and is plugged into a grounded socket.
- **Keep Away from Children:** The device and its accessories should be stored securely out of reach.
- **Report Issues Promptly:** If you notice abnormal functioning, contact your distributor or an authorized technician immediately.

Safety Instructions

- Plug the device into a properly grounded electrical outlet with the appropriate voltage.
- Do not use the device if it or its power cable is damaged.
- Never touch the device, power plug, or cables with wet hands.
- Avoid attempting to repair or modify the device. Repairs should only be performed by authorized technicians.
- Store the device in a cool, dry place away from direct sunlight and moisture.




Warnings (High Frequency Mode)

- When using the high-frequency electrotherapy, customers must not touch any metal objects, any and all jewelry or electronics must be removed.
- When using high-frequency electrotherapy, do not use skin care creams or masks containing alcohol. If you must use a toner containing alcohol, apply it after using high-frequency electrotherapy.
- **Never** point any type of electrotherapy tube directly at the eyeball.
- Do not touch the high-frequency electrotherapy device output interface or handle output port on the main unit.
- This device should only be operated by trained beauty professionals. The attachments generate high-frequency electricity, so keep them away from other electrical devices and use caution.


Warnings (Galvanic Mode)

- Do not allow the export and import probes to remain in contact with each other for long periods of time.
- The galvanic device requires ample lubrication during use; otherwise, the effects of the treatment will not be done properly.
- If the roller probe stops rotating, please refrain from using it and replace the attachment.
- When using the device, do not splash water or other corrosive liquids onto the main unit.
- Solution Compatibility: Match the skincare solution to the selected polarity. When the probe is positive (+), use an acidic solution; when the probe is negative (–), use an alkaline solution. If the solution type is unknown, use pH test strips to determine compatibility.
- It is strictly prohibited to use ball probes and roller probes with damaged, rusted, or scratched surfaces, as this may cause skin abrasions. Probes must be carefully protected and must not be scratched or damaged by other hard objects.

Included Items

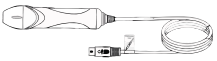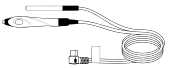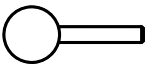Qty	Part Name	Image
1	High Frequency + Galvanic Unit	
1	Accessory Rack	
1	High Frequency Electrotherapy Wand	
1	Mushroom Electrode	
1	Spoon-Shaped Electrode	
1	Dot Electrode	
1	Straight Electrode	
1	Galvanic Handle	
1	Roller Probe	
1	Round Eye Probe	
1	Galvanic Sphere	
1	Power Cord	

High Frequency Tool Components

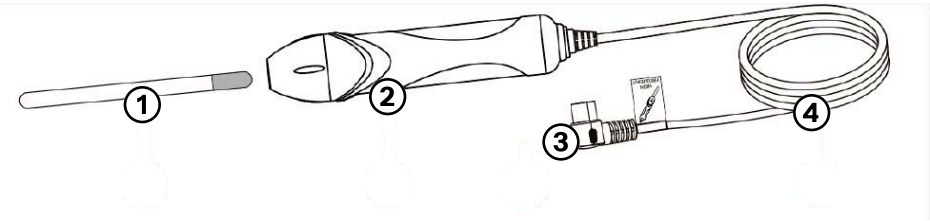


1. Electrotherapy Tube
2. Electrotherapy Handle
3. Plug for Handle
4. Flexible Cable

Selecting The Right Glass Electrotherapy Tube

- **Mushroom Shape:** This tube is suitable for treating large areas such as the cheeks, forehead, and neck.
- **Spoon Shaped:** Ideal for medium-sized areas such as the chin.
- **Dot Shape:** Designed for treating small areas like the nasal pit.

Galvanic Tool Components

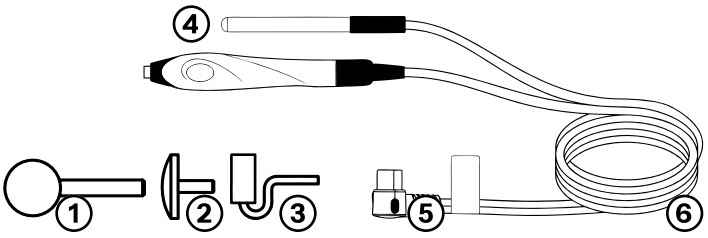

1. Round Eye Probe
2. Galvanic Sphere
3. Roller Probe
4. Conduction Probe
5. Conduction Probe Plug
6. Conduction Probe Cable

Selecting The Right Galvanic Probe:


- **Round Probe:** Best for around the eye area.
- **Sphere Probe:** Best for face, everywhere except the eye area.
- **Roller Probe:** Best for body and face, covers larger areas.
- **Conduction Probe:** The device and its accessories should be stored securely out of reach.

Touch Panel Interface




1. Time display
2. Decrease Timer
3. Increase Timer
4. Galvanic: Normal
5. Galvanic: Sensitive
6. Galvanic: Extra Sensitive
7. Galvanic Function Toggle
8. Play/Pause Button
9. High Frequency Toggle
10. Anion Introduction
11. Cation Introduction
12. Decrease Power Level
13. Increase Power Level
14. Power Level Display

Device Components

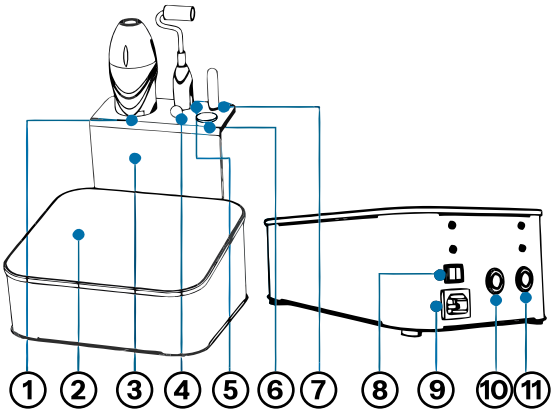


1. High Frequency Handle
2. Touch Panel
3. Accessory Rack
4. Galvanic: Round Probe
5. Galvanic: Roller Probe
6. Galvanic: Sphere Probe
7. Galvanic: Grip
8. Power Switch
9. AC Power Cable Socket
10. Galvanic Socket
11. High Frequency Socket

Assembly Instructions

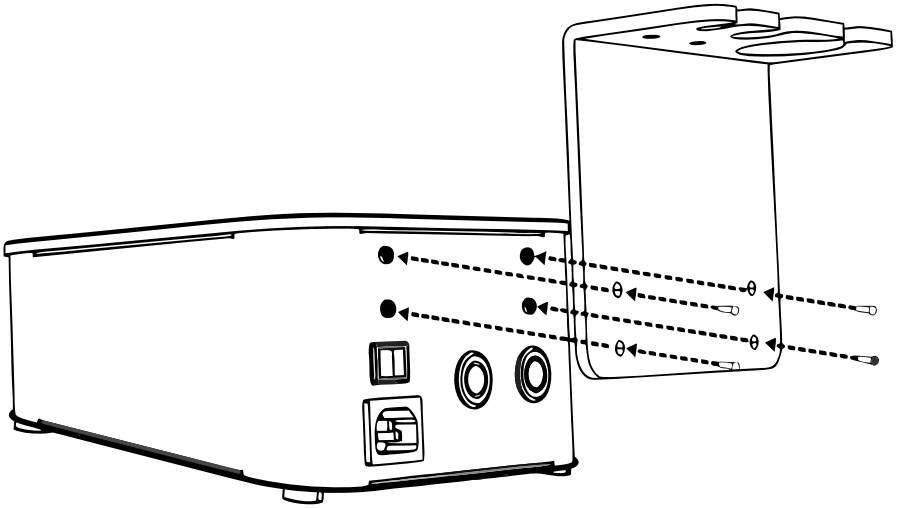

- **Attach the accessory rack:** Remove the pre-attached screws and line up the acrylic rack and secure it to the unit with the screws.
- **Place the probes and wand into the rack.**
- **Plug the cables into their corresponding ports.**

Operation Instructions

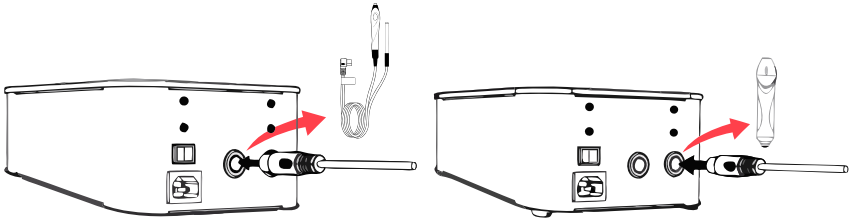

- **Flip the power switch on to enter standby mode:** The time display will show “15,” and the energy display will show “1.”
- In standby mode, press the Start/Pause button to begin operation. Settings can be adjusted while in standby mode.
- Press the Timer button, then use the + and – buttons to adjust the operating time. The adjustable range is 1 to 30 minutes.
- Press the Level button, then use the + and – buttons to adjust the energy output. The adjustable range is Level 1 to Level 9.

High Frequency Information




High-frequency energy ionizes the gas inside the glass electrodes, creating ultraviolet rays and a glowing effect. This current stimulates the skin and provides therapeutic benefits.

- **Uses Include:**
  - Stimulate glandular activity.
  - Enhance blood circulation.
  - Provide germicidal and antibacterial effects.
  - Promote deep product penetration.
  - Gently warm the skin for relaxation and enhanced absorption of skincare products.
  - This device is ideal for professional estheticians seeking to offer high-quality treatments with proven results.

High Frequency Operating Instructions

- **Pre-Treatment Prep:**
  a. Remove all jewelry from both the esthetician and the client.
  b. Ensure the client is not in contact with metal objects (e.g., the treatment chair frame).
  c. Have the client remove contact lenses before treatment.
  d. Explain the procedure to the client, including the sparking sensation they may feel.

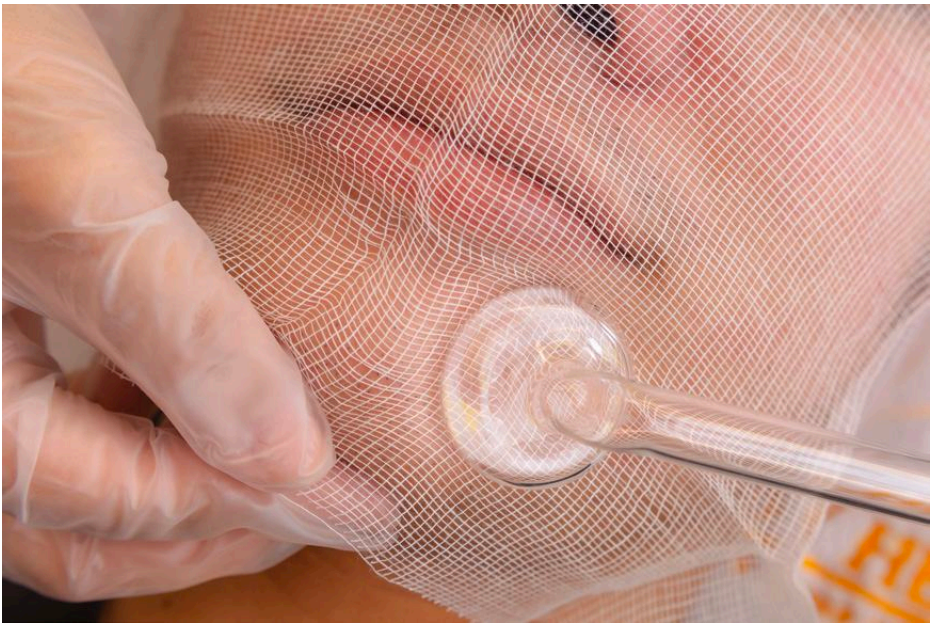

- **Direct Technique (Soothing):**
  a. Apply a non-flammable, water-soluble product or gauze to the skin for slip. Avoid alcohol-based products.
  b. Set the current intensity to its lowest setting before turning the unit “On” by pressing the play button.
  c. Insert the chosen electrode into the handle, ensuring it is free of dust or debris.
  d. Place your finger on the electrode, turn the unit “On,” and position the electrode on the client’s skin. Remove your finger from the electrode.
  e. Adjust the intensity based on the client’s comfort level.
  f. Use small, circular movements starting from the neck and moving upward to the jaw, cheeks, nose, and forehead.



- **Indirect Technique (Soothing):**
  a. Have the client hold the straight electrode securely in both hands.
  b. Cover the electrode with a gauze or cheesecloth for comfort.
  c. Ensure the client maintains a firm grip before turning the unit “On” with the play button.
  d. Place one hand on the client’s skin at the treatment area, then proceed with massage techniques.




Galvanic Operating Instructions

1. Turn on the main power switch on the back of the device. The power light on the front will turn on, and the time display will show “15”. The energy display window will display “1”.
2. Select Polarity (Manual Control): Before starting treatment, manually select the required ion polarity based on your product's pH:
   - Press the Cation (+) for positive ion introduction.
   - Press the Anion (-) for negative ion introduction or deep cleansing.
3. Use the Level Increase/Decrease buttons to adjust the energy output. The adjustable range is Level 1 to Level 9.
4. Once the probe is in contact with the client, press the Start/ Pause button to begin flow of current.

Ion Introduction Polarity Selection

- **(+): Cation (positive ion) Introduction**
- **(-): Anion (negative ion) Introduction**

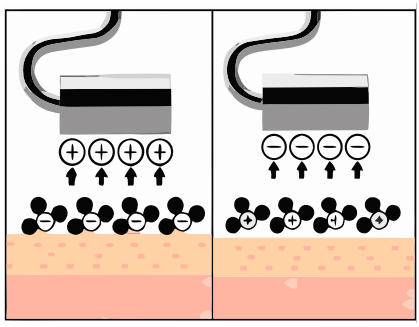


Usage (Iontophoresis Function)

A. Monopolar Infusion Method

1. Apply half of the nutrient solution directly to the skin. Soak the remaining solution into a cotton pad and wrap it around the probe.
2. Set the Level to minimum. Have the client hold the metal electrode rod while the operator holds the probe and places it on the treatment area.
3. Press the Start/Pause button to activate the device.
4. Select the required ion polarity (+ or −) based on the product used.
5. Gradually increase the Level until a mild tingling sensation is felt.
6. Move the probe slowly in a zigzag or circular motion, maintaining constant contact with the skin. (Figure 1)
   - Typical treatment time: 3–5 minutes

Figure 1

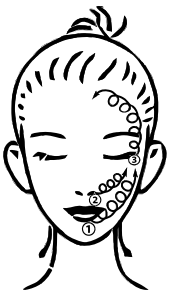

7. After completion:
   - Press Start/Pause to stop output
   - Turn off the Power
   - Remove the probe
   - Have the client release the electrode rod

B. Impurity Export (Deep Cleansing)

1. Set the working probe to negative (−) polarity. The client holds the positive (+) electrode rod..
2. Soak a cotton pad with saline or impurity-removal solution and wrap it around the probe, or apply gauze/mask to the treatment area.
3. Turn on the device and select anion (−) mode.
4. Start at a low energy level and increase gradually.
5. Move the probe slowly across the treatment area, keeping full contact with the skin.
6. After treatment, reduce energy to minimum, stop the device, and remove the probe.

C. Alternating Output Mode

1. Follow the same procedure as Method A (Monopolar Infusion).
2. In this mode, the device automatically alternates between positive and negative currents.
3. The polarity switches approximately every 0.5 seconds, providing continuous electrical stimulation to the skin. The rapid switching between positive and negative charges creates a "micro-massage" effect at the cellular level, which can help wake up sluggish or tired skin.

Specifications

Product Name:	2 in 1 High Frequency Galvanic Elite Edition
Voltage:	AC100-120V/220-240V
Frequency:	50/60Hz
Power Consumption:	20W
Dimensions(L*W*H):	12x10x9 in
Net Weight:	4.92 lbs

Warranty Information

The manufacturer (the “warrantor”) warrants to the original purchaser that this equipment will be free from defects under normal and proper operation for a period of one (1) year from the date of purchase. This warranty covers all parts and labor for that one (1) year period. This warranty shall not apply to any equipment that has been subjected to abuse, misuse and alteration of any type or to any defect or damage caused by repair, replacement, substitution, or use with parts other than parts approved by the manufacturer. The manufacturer shall not be liable for the loss of use of this equipment, loss of time or business, inconvenience or any other indirect, consequential, special or incidental damages due to defects covered by this warranty or any implied warranty. To implement this warranty, contact the dealer from whom you purchased the equipment. If necessary, you will be asked to return all, or the appropriate part, of the equipment. If warranted, the repairs will be performed at no charge and the equipment returned to you. Unit owner must prepay all shipping costs when returning items for warranty service. This limited warranty is the only written or express warranty given by the manufacturer. This warranty gives you specific legal rights but you may also have other rights, which may vary from state to state or province to province.

SKIN+ACT

Phone: 1-800-418-0830
Email: info@skinact.com
Website: skinact.com